



Exo-Check Exosome Antibody Arrays

Cat # EXORAY200A-4, Cat #EXORAY210A-8

User Manual

Please see individual components for storage conditions

Version 3
8/21/2018

A limited-use label license covers this product. By use of this product, you accept the terms and conditions outlined in the License and Warranty Statement contained in this user manual.

Contents

Product Description.....	2
List of Components.....	2
Storage.....	2
General Information	2
Protocol for Exo-Check Arrays	3
Next Steps and Related Products	4
Example Data and Applications	4
Technical Support	6
Licensing and Warranty Statement	7

Product Description

The Exo-Check antibody array (4 and 8 array formats) has 12 pre-printed spots including 8 antibodies for known exosome markers (CD63, CD81, ALIX, FLOT1, ICAM1, EpCam, ANXA5 and TSG101) and 4 controls. The two positive control signals provide a bright signal that indicates that the detection reagents are working properly. The Blank spot serves as a background control. The GM130 cis-Golgi marker monitors any cellular contamination in your exosome isolations. The array uses protein lysate from exosomes obtained from a wide variety of methods, including ultracentrifugation, precipitation-based, and size exclusion chromatography. The kits come complete with a secondary detection mixture conjugated to HRP. The Exo-Check arrays are semi-quantitative and can be used to evaluate the relative abundance of certain exosome protein markers from a given set of samples.

List of Components

Components	4 array kit	8 array kit	Storage Temperature
Lysis Buffer (10X)	1 ml	2 ml	4°C
Labeling Reagent	4 x 20 µl	8 x 20 µl	-20°C
Columns	4	8	4°C
Column Buffer	20 ml	40 ml	4°C
Blocking Buffer	20 ml	40 ml	4°C
Wash Buffer (20X)	5 ml	10 ml	4°C
Detection Buffer	20 ml	40 ml	4°C
Membranes	4	8	4°C

Not supplied: HRP developer reagent. SBI recommends the use of **Advansta WesternBright Sirius HRP substrate, catalog# K12043-C20.**

Storage

The kits are shipped on blue ice and should be stored at +4°C (**except for Labeling Reagent, which must be stored at -20°C**) upon receipt. Properly stored kits are stable for 3 months from the date received.

General Information

The various exosome antibody spots will provide signals of varying degree, depending upon the source of the isolated exosomes. Note that for a given sample, signals from each spot on the array cannot be directly compared

since the antibody concentrations are different across the array. However, it can be compared across samples to provide a qualitative assessment of a given spot's signal intensity.

Protocol for Exo-Check Arrays

Sample preparation

1. Start with isolated EVs resuspended in 1X PBS.
2. Perform a protein quantitation assay (e.g. BCA or Qubit assay) to determine the amount of protein in your sample. Transfer 50 µg of sample to a new tube; save the rest for future use.
3. Add 10X Lysis Buffer to the 50 µg of sample (from step 2) to a final concentration of 1X. Vortex for 15-30 seconds.

Note: No centrifugation is needed during the lysate preparation.

4. Equilibrate the Labeling Reagent to room temperature.
5. Add 1 µl of Labeling Reagent* to the entire content from step 3. Vortex to mix the sample.
*Use one amber tube of Labeling Reagent per reaction/membrane. Do not save the reagent. Always use a fresh amber tube for preparing a new sample.
6. Incubate at room temperature for 30 mins with constant mixing.
7. Remove excess Labeling Reagent using the provided columns:
 - a. Vortex the column to resuspend the medium.
 - b. Loosen the screw cap and twist off the bottom closure.
 - c. Place the column in a collection tube.
 - d. Remove the storage buffer by centrifuging at 800 x g for 1 min. Discard the storage buffer.
 - e. Equilibrate the column by adding 400 µl Column Buffer.
 - f. Centrifuge for 1 min at 800 x g. Discard the flow-through and place the column into the collection tube.
 - g. Repeat steps e-f 4 times for a total of 5 times.
 - h. Apply the sample from step 6 (100-180 µl) slowly in the middle of the packed bed.
Note: If the sample was less than 140 µl, add a stacker volume of Column Buffer to reach 140 µl.
 - i. Elute the sample (i.e., labeled exosome lysate) by centrifuging at 800 x g for 2 mins.
8. In a new 15 ml conical tube, combine the 50 µg of labeled exosome lysate with 5 ml Blocking Buffer and mix by inverting the tube 3 times.

Antibody array capture

1. Remove the membrane from the white envelope using forceps, grabbing the edge of the membrane.
2. In a small, clean tray, briefly wet the membrane in 5 ml distilled water for 2 mins at room temperature.
3. Carefully, decant the water from the membrane. Pipette the 5 ml labeled exosome lysate/blocking buffer mixture to one membrane. Place the membrane "face-up" by positioning the membrane such that the upper left notched corner can be seen in the left-hand side on top.
4. Incubate the tray/membrane mixture shaking at 2-8°C overnight on a shaker or rocker.

Membrane washing and signal detection

1. Next day, decant the lysate/blocking mixture from the membrane carefully.
2. Add 5 ml of 1X Wash Buffer (diluted from 20X stock with DiH2O) and rock gently for 5 min at room temperature.
3. Decant 1X Wash Buffer. Repeat step 2 once more. Decant after the last wash.
4. Add 5 ml Detection Buffer to the membrane and incubate for 30 min at room temperature on a shaker or rocker.
5. Prepare the Developer mixture (not included) about 5 mins before the end of step 4 according to the manufacturer's recommendation.
6. Decant the Detection Buffer and add 5 ml 1X Wash Buffer and shake for 5 min at room temperature.
7. Decant the 1X Wash Buffer, repeat wash twice for a total of three washes. Decant after the last wash.
8. Develop and image the membrane according to instrument manufacturer's recommendations.

Next Steps and Related Products

Application	Related Products	Website links
Precipitation of Exosomes from other biological fluids		
Exosome Isolation from Tissue Culture Media	ExoQuickTC	https://www.systembio.com/microrna-research/exoquick-exosomes/ordering
Exosome Isolation from Plasma	ExoQuick Plasma prep and Exosome precipitation kit	https://www.systembio.com/microrna-research/exoquick-exosomes/ordering
Protein Characterization of Exosomes		
Western blotting	Exosome antibodies	https://www.systembio.com/microrna-research/exosome-antibody/exoab
ELISA	ExoELISA Kits	https://www.systembio.com/microrna-research/exosome-antibody/elisas
Quantification of Exosomes		
Quantification of exosomes	EXOCET Assays	https://www.systembio.com/microrna-research/exosome-antibody/exocet-assay
RNA extraction from Exosomes		
RNA extraction and profiling	SeraMir kits	https://www.systembio.com/microrna-research/seramir-exosome-rna-profiling/overview

Example Data and Applications

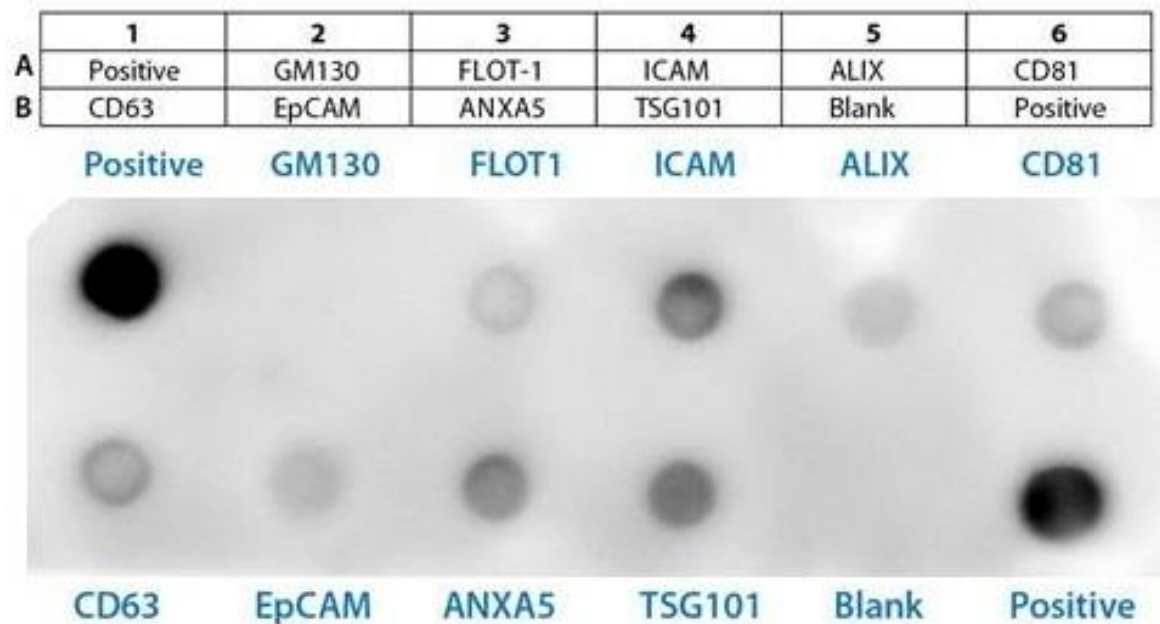


Figure 1. Data for array spots obtained from 50 μ g of normal pooled human serum exosome (isolated using ExoQuick) lysate.

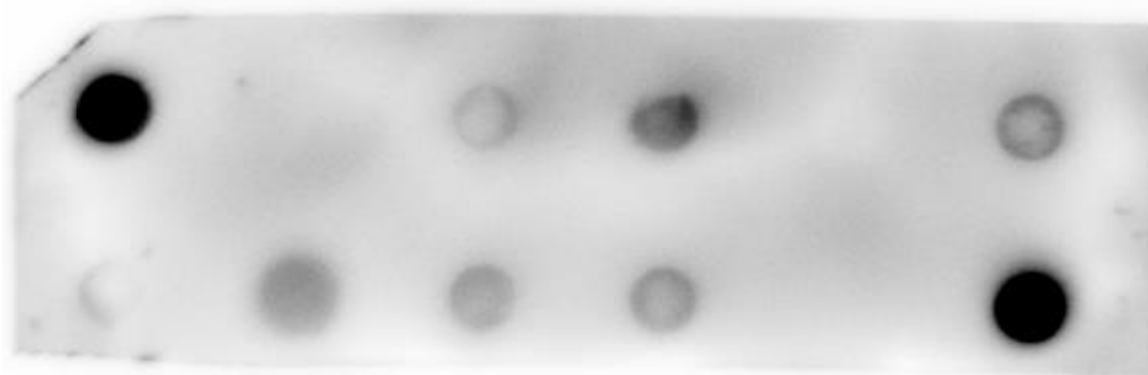


Figure 2. Data for array spots obtained from 100 μ g of HEK293T exosome (isolated using ExoQuick-TC) lysate. Key for individual spots is same as in Figure 1.

Information about the exosome protein targets

SPOT POSITION*	ID	NOTES	LINK TO EXOCARTA ENTRY FOR THE HUMAN PROTEIN
A1	Positive control	Positive control for HRP detection (derived from human serum exosomes)	NA
A2	GM130	Cis-golgi matrix protein—control for cellular contamination in exosome preparation	NA
A3	FLOT1	Flotillin-1	ExoCarta_10211
A4	ICAM1	Intercellular adhesion molecule 1	ExoCarta_3383
A5	ALIX	Programmed cell death 6 interacting protein (PDCD6IP)	ExoCarta_10015
A6	CD81	Tetraspanin	ExoCarta_93185
B1	CD63	Tetraspanin	ExoCarta_967
B2	EpCam	Epithelial cell adhesion molecule; often found in cancer-derived exosomes	ExoCarta_4072
B3	ANXA5	Annexin A5	ExoCarta_308
B4	TSG101	Tumor susceptibility gene 101	ExoCarta_7251
B5	Blank	Background control	NA
B6	Positive control	Positive control for HRP detection (derived from human serum exosomes)	NA

Table 1. Information for each spot on the array

Technical Support

For more information about SBI products and to download manuals in PDF format, please visit our web site: <http://www.systembio.com>

For additional information or technical assistance, please call or email us at:

System Biosciences (SBI)
2438 Embarcadero Way
Palo Alto, CA 94303

Phone: (650) 968-2200
(888) 266-5066 (Toll Free)

Fax: (650) 968-2277

E-mail:

General Information: info@systembio.com
Technical Support: tech@systembio.com
Ordering Information: orders@systembio.com

Licensing and Warranty Statement

Limited Use License

Use of the Exo-Check Exosome Antibody Arrays (*i.e.*, the “Product”) is subject to the following terms and conditions. If the terms and conditions are not acceptable, return all components of the Product to System Biosciences (SBI) within 7 calendar days. Purchase and use of any part of the Product constitutes acceptance of the above terms.

The purchaser of the Product is granted a limited license to use the Product under the following terms and conditions:

- The Product shall be used by the purchaser for internal research purposes only. The Product is expressly not designed, intended, or warranted for use in humans or for therapeutic or diagnostic use.
- The Product may not be resold, modified for resale, or used to manufacture commercial products without prior written consent of SBI.
- This Product should be used in accordance with the NIH guidelines developed for recombinant DNA and genetic research.

SBI has pending patent applications related to the Product. For information concerning licenses for commercial use, contact SBI.

Purchase of the product does not grant any rights or license for use other than those explicitly listed in this Licensing and Warranty Statement. Use of the Product for any use other than described expressly herein may be covered by patents or subject to rights other than those mentioned. SBI disclaims any and all responsibility for injury or damage which may be caused by the failure of the buyer or any other person to use the Product in accordance with the terms and conditions outlined herein.

Limited Warranty

SBI warrants that the Product meets the specifications described in this manual. If it is proven to the satisfaction of SBI that the Product fails to meet these specifications, SBI will replace the Product or provide the purchaser with a refund. This limited warranty shall not extend to anyone other than the original purchaser of the Product. Notice of nonconforming products must be made to SBI within 30 days of receipt of the Product.

SBI’s liability is expressly limited to replacement of Product or a refund limited to the actual purchase price. SBI’s liability does not extend to any damages arising from use or improper use of the Product, or losses associated with the use of additional materials or reagents. This limited warranty is the sole and exclusive warranty. SBI does not provide any other warranties of any kind, expressed or implied, including the merchantability or fitness of the Product for a particular purpose.

SBI is committed to providing our customers with high-quality products. If you should have any questions or concerns about any SBI products, please contact us at (888) 266-5066.

This page intentionally left blank

This page intentionally left blank

This page intentionally left blank



System Biosciences (SBI)
2438 Embarcadero Way
Palo Alto, CA 94303

Phone: (650) 968-2200
Toll Free (888) 266-5066
Fax: (650) 968-2277

E-mail:

General Information: info@systembio.com
Technical Support: tech@systembio.com
Ordering Information: orders@systembio.com