



---

# EVery cDNA Synthesis Kit

Cat # EVery200B-1

## User Manual

---

**Storage:**

Store at -20°C upon receipt

Version 1  
10/10/2022

A limited-use label license covers this product. By use of this product, you accept the terms and conditions outlined in the License and Warranty Statement contained in this user manual.



---

## Contents

Product Description .....	2
List of Components .....	3
Additional Required and Optional Equipment Not Included in Kit .....	3
Protocol .....	4
Example Data and Applications .....	5
Technical Support .....	8
Licensing and Warranty Statement .....	8



# Product Description

Get your EV RNA ready for downstream analysis with this quick and easy cDNA synthesis kit optimized for use with the [EVery EV RNA Isolation Kit](#).

- Get quick (<2 hrs) and easy cDNA synthesis optimized for use with EVery EV RNA Isolation Kit
- Use with most downstream applications such as PCR, qPCR and miRNA profiling
- Includes sufficient Universal reverse primer for all 8 profiles from EVery miRNome Profiler for Human Serum and Plasma, EVery500B-1

Optimized for use with [EVery EV RNA Isolation Kit](#), EVery cDNA Synthesis Kit provides high quality exosomal cDNA for your downstream applications such as PCR, qPCR, miRNA profiling, and RNA-seq. The kit generates cDNA through poly(A) tailing of exosomal RNA followed by reverse transcription with MMLV Reverse Transcriptase. This method uses the template-switching oligonucleotide property of MMLV RT to generate cDNA molecules with two adaptor primers from both 5'- and 3'- ends (Figure 1), the expression of which can be detected by qPCR using EVery miRNome Universal Reverse Primer included in the kit. EVery cDNA Synthesis Kit comes with sufficient reagents to perform 20 reactions.

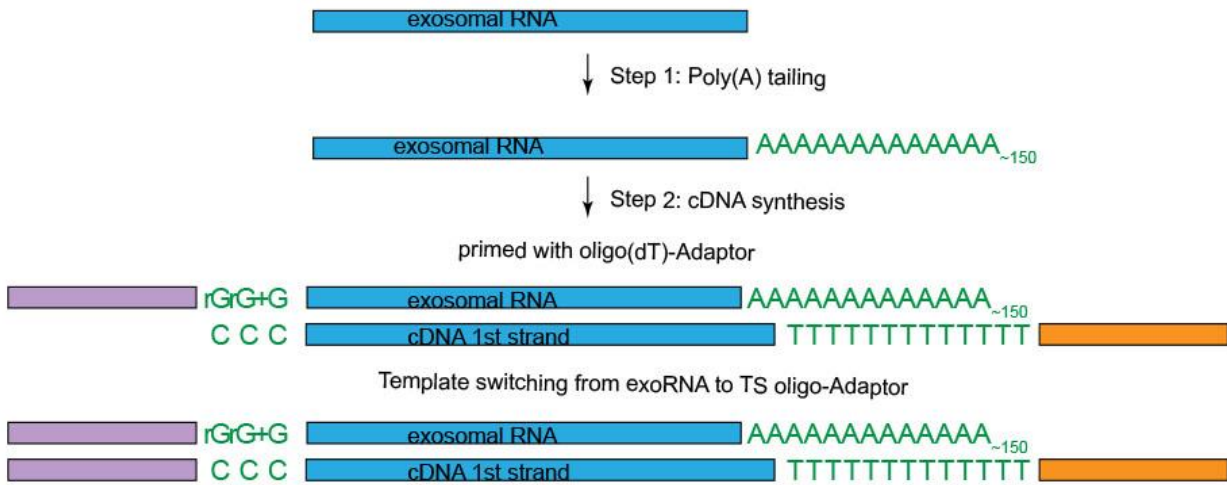


Figure 1. EVery cDNA synthesis workflow.

## List of Components

**Table 2. Components of EVery cDNA Synthesis Kit, EVery200B-1**

Components	Qty/Volume	Storage Temperature
10X Poly(A) tailing reaction buffer	20 µL	-20°C
10 mM ATP	20 µL	-20°C
Poly(A) tailing enzyme	10 µL	-20°C
RNase inhibitor	20 µL	-20°C
RNase free water	400 µL	-20°C
Oligo(dT)-Adaptor	40 µL	-20°C
5X First strand buffer	100 µL	-20°C
DTT	40 µL	-20°C
dNTPs	40 µl	-20°C
TS oligo-Adaptor	20 µL	-20°C
First strand synthesis enzyme	40 µL	-20°C
Universal Reverse Primer	3 vials of 500 µL	-20°C

### Additional Required and Optional Equipment Not Included in Kit

1. PCR tubes
2. Thermal cycler
3. EVery miRNA Spike-In Kit, EVery600B-1 can be used to determine the efficiency of cDNA synthesis.

## Protocol

### Poly(A) tailing reaction and cDNA synthesis:

1. Set up the following reaction:

Component	Quantity
10X Poly(A) tailing reaction buffer	1 $\mu$ L
10 mM ATP	1 $\mu$ L
Poly(A) tailing enzyme	0.5 $\mu$ L
RNase inhibitor	0.25 $\mu$ L
EVs RNA	X $\mu$ L
RNase free water	Up to 10 $\mu$ L

2. Gently mix the reaction and centrifuge briefly.
3. Incubate the tube at 37°C for 20 minutes.
4. Add **2  $\mu$ L of Oligo(dT)-Adaptor** to prime the cDNA.
5. Incubate at 72°C for 3 minutes and then, 2 minutes on ice.
6. Add the following components for the first strand synthesis reaction:

Component	Quantity
5X First strand buffer	5 $\mu$ L
dNTPs	2 $\mu$ L
DTT	2 $\mu$ L
TS oligo- Adaptor	1 $\mu$ L
First strand synthesis enzyme	2 $\mu$ L
RNase free water	Up to 25 $\mu$ L

7. Place the tube in the thermal cycler with a heated lid, preheated to 42°C. Run the following program:

42°C – 90 min

70°C – 10min

4°C – forever

## Example Data and Applications

### Get robust cDNA synthesis with EVery cDNA Synthesis Kit

To demonstrate the excellent RNA yields and robust cDNA synthesis obtained with EVery family of products, we isolated EVs from 250  $\mu$ L of serum using [SmartSEC Single](#), spiked in 0.1 pmol of Cel-miR-39, and used both [EVery EV RNA Isolation Kit](#) and a phenol-based kit to isolate RNA. The isolated RNA was reverse transcribed using EVery cDNA Synthesis Kit and the copy number of Cel-miR-39 measured (Figure 2). EVery EV RNA Isolation Kit delivered similar levels of Cel-miR-39 as the phenol-based method.

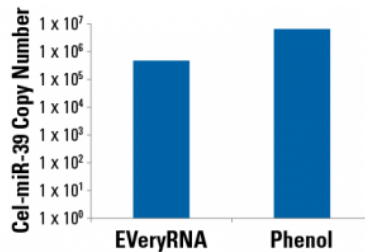


Figure 2. EVery is every bit as good as phenol.

### EVERy efficiently isolates mRNA

We used EVery EV RNA Isolation Kit and EVery cDNA Synthesis Kit to isolate mRNA and synthesize cDNA from cells overexpressing eGFP (Figure 3). Robust levels of eGFP mRNA are recovered and converted to cDNA when cells are overexpressing eGFP.


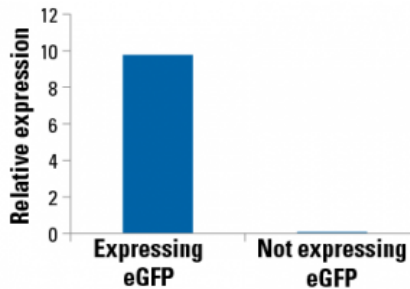


Figure 3. EVery efficiently isolates longer RNAs like mRNA




**miRNA isolated with EVery100B-1 and converted to cDNA using EVery200B-1 can be used for miRNA profiling**

We isolated EVs from 250 µL of normal or breast cancer patient serum using [SmartSEC Single](#) and used EVery EV RNA Isolation Kit and EVery cDNA Synthesis Kit to isolate and reverse transcribe EV RNAs for miRNA profiling using [EVery miRNome Profiler for Human Serum and Plasma \(EVery500B-1\)](#). Distinct miRNA expression profiles of breast cancer patient serum compared to normal serum are shown in Figure 4. Differential expression of selected breast cancer associated miRNA markers were further confirmed by qPCR (Figure 5).

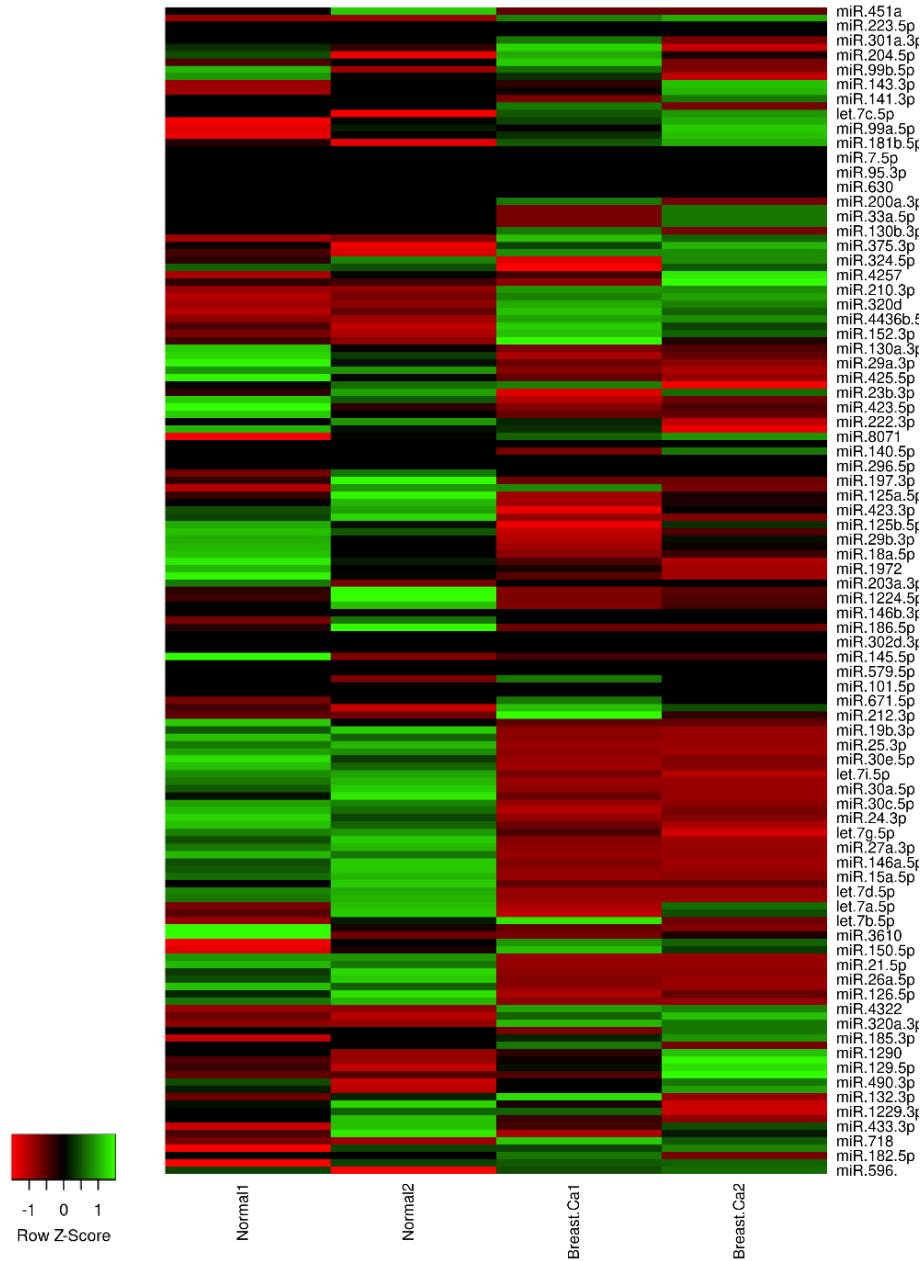

**Figure 4. EVery family of products generate high-quality cDNA suitable for miRNA profiling**

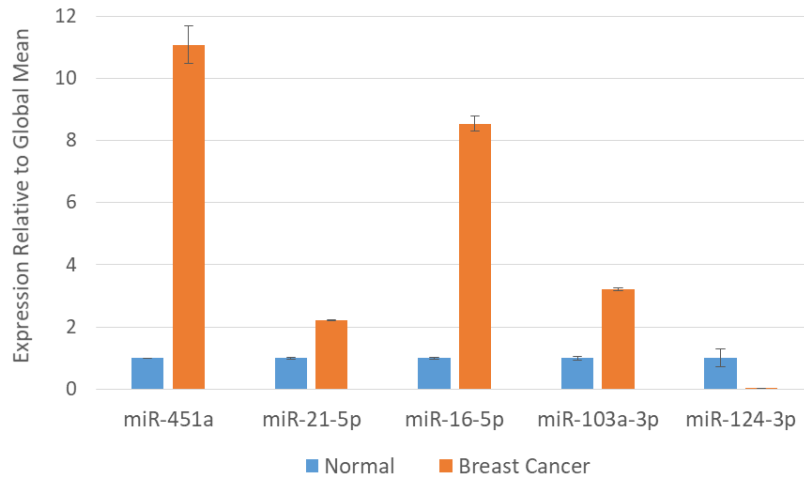

**Figure 5. Differential expression of selected miRNAs in normal vs. breast cancer serum confirmed by qPCR. Data normalized with global mean. Expression of miR-451a, miR-21-5p, miR-16-5p and miR-103a are significantly upregulated, while expression of the tumor suppressor miR-124-3p is significantly downregulated in breast cancer serum compared to normal serum.**

## Technical Support

For more information about SBI products and to download manuals in PDF format, please visit our web site: <http://www.systembio.com>

For additional information or technical assistance, please call or email us at:

System Biosciences (SBI)  
2438 Embarcadero Way  
Palo Alto, CA 94303

**Phone:** (650) 968-2200  
**Toll-Free:** (888) 266-5066  
**Fax** (650) 968-2277

**E-mail:**  
General Information: [info@systembio.com](mailto:info@systembio.com)  
Technical Support: [tech@systembio.com](mailto:tech@systembio.com)  
Ordering Information: [orders@systembio.com](mailto:orders@systembio.com)

## Licensing and Warranty Statement

### Limited Use License

Use of EVERY EV RNA Isolation Kit (*i.e.*, the “Product”) is subject to the following terms and conditions. If the terms and conditions are not acceptable, return all components of the Product to System Biosciences (SBI) within 7 calendar days. Purchase and use of any part of the Product constitutes acceptance of the above terms.

The purchaser of the Product is granted a limited license to use the Product under the following terms and conditions:

- The Product shall be used by the purchaser for internal research purposes only. The Product is expressly not designed, intended, or warranted for use in humans or for therapeutic or diagnostic use.
- The Product may not be resold, modified for resale, or used to manufacture commercial products without prior written consent of SBI.
- This Product should be used in accordance with the NIH guidelines developed for recombinant DNA and genetic research.

This Product is covered by licensed Patents. For information concerning licenses for commercial use, contact SBI.

Purchase of the product does not grant any rights or license for use other than those explicitly listed in this Licensing and Warranty Statement. Use of the Product for any use other than described expressly herein may be covered by patents or subject to rights other than those mentioned. SBI disclaims any and all responsibility for injury or damage which may be caused by the failure of the buyer or any other person to use the Product in accordance with the terms and conditions outlined herein.

### Limited Warranty

SBI warrants that the Product meets the specifications described in this manual. If it is proven to the satisfaction of SBI that the Product fails to meet these specifications, SBI will replace the Product or provide the purchaser with a refund. This limited warranty shall not extend to anyone other than the original purchaser of the Product. Notice of nonconforming products must be made to SBI within 30 days of receipt of the Product.

SBI's liability is expressly limited to replacement of Product or a refund limited to the actual purchase price. SBI's liability does not extend to any damages arising from use or improper use of the Product, or losses associated with the use of additional materials or reagents. This limited warranty is the sole and exclusive warranty. SBI does not provide any other warranties of any kind, expressed or implied, including the merchantability or fitness of the Product for a particular purpose.

SBI is committed to providing our customers with high-quality products. If you should have any questions or concerns about any SBI products, please contact us at (888) 266-5066.

© 2020 System Biosciences (SBI), All Rights Reserved

***This page intentionally left blank.***

*This page intentionally left blank.*

***This page intentionally left blank.***



---

System Biosciences (SBI)  
2438 Embarcadero Way  
Palo Alto, CA 94303

**Phone:** (650) 968-2200  
**Toll-Free:** (888) 266-5066  
**Fax:** (650) 968-2277

**E-mail:**

General Information: [info@systembio.com](mailto:info@systembio.com)  
Technical Support: [tech@systembio.com](mailto:tech@systembio.com)  
Ordering Information: [orders@systembio.com](mailto:orders@systembio.com)